



Bariatric Wheelchair Minimaxx With Push Motor

Strongly engineered and functionally designed exclusively for bariatric users.

Bariatric Wheelchair Minimaxx With Push Motor is a strong power pack that assists the attendant pushing the wheelchair with no effort. Using a Minimaxx With Push Motor the caregiver avoids sore shoulders and stress on wrists.

Even on slopping surface up to 5% Minimaxx With Push Motor proudly helps over you over the top and going down you find the security of our disc braking system. If an immediate stop is needed just pull the lever and the wheelchair stops at once.

The push motor is easy to operate. The control panel is placed on the handle and offers Drive and Reverse mode as well as speed control from 1-5 km/h.

As for the construction of the wheelchair STW (Stiffness To Weight) was our goal. We used very light tubing technology as used on high end MTB and road bikes. That is why Minimaxx With Push Motor only weighs approx. 9% of its weight capacity and still provides a smooth drive thanks to the light tubing.

The castor wheels are further in front as on standard wheelchairs. This prevents the wheelchair from tipping.

The very long armrests make it easy for the user to get in and out of Minimaxx. The low seat height makes it easy to access the wheelchair.

The footrests are designed to carry 100 kg each. The footrests are also side and width adjustable. This enables users with legs aligned or with a wide leg spread to use the wheelchair.

MAX 325 kg

- Push motor is easy to operate
- Different armrests available



Push motor



Standard armrests



Sport armrests




Sport armrests, side view

Bariatric Wheelchair Minimaxx With Push Motor Specifications

Wheelchair Minimaxx with Push Motor	
Maximum weight capacity is	325 kg
Seat width	56 / 61 / 66 / 71 cm
Total width excluding push handle	80 / 85 / 90 / 95 cm
Total height including push handle, lowest position	110 cm
Total height	103 cm
Total length without footrests	70 cm
Seat depth	46 cm
Seat height	43 cm
Backrest height	51 cm
Weight	Approx. 27 kg
Weight push motor	23 kg
Weight battery	9 kg
Battery	12V 12AH x 2
Motor	24V 450W 4.5 amps
Tire size	50 mm x 200 mm (Solid Tire) x 2
Charger capacity	Off Board (1.8A)
Range per charge	4.5 km
Max. speed	5 km/h
Gradient	5°
ISO - classification	ISO 12 23 12 - 02

Ideal use for...

- Bariatric users in nursing homes, hospitals or private homes.

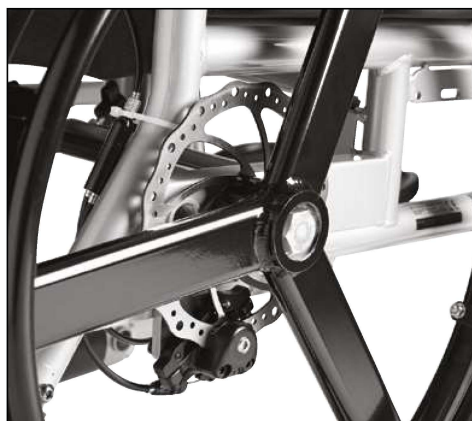

Disc brakes

Model	Seat width cm	Seat depth cm	Order No
XXL-Rehab Minimaxx With Push Motor, incl. disc brakes and footrests	56	46	0110-056-000
XXL-Rehab Minimaxx With Push Motor, incl. disc brakes and footrests	61	46	0110-061-000
XXL-Rehab Minimaxx With Push Motor, incl. disc brakes and footrests	66	46	0110-066-000
XXL-Rehab Minimaxx With Push Motor, incl. disc brakes and footrests	71	46	0110-071-000

Accessories	Order No
Sport armrests, pair	0100-999-006
Table	0100-999-001
IV holder and O2 holder, diameter 10 cm	0100-999-020
Footrest swing-away, right	0100-999-010
Footrest swing-away, left	0100-999-011



Table



Hand brake operated by attendant



IV holder and O2 holder



Sport armrests