




TIGER OUTDOOR

FOREARM WALKER

TIGER OUTDOOR FOREARM WALKER

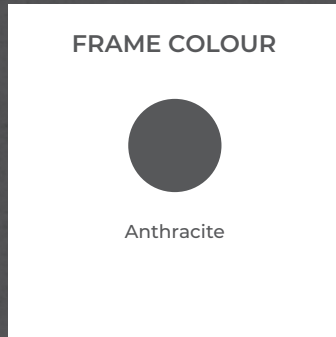

The Mobilex outdoor Tiger forearm walker is based on the Lion off road rollator and as such is designed for use anywhere. It is perfect for users that require the additional support of the forearm rests, either as a result of lack of wrist strength/support or rehabilitation purposes. Although very robust, it is still relatively lightweight at just 12kgs and offers a high degree of stability and ease of movement.

The forearm rests are height-adjustable, with an easy to read numerical marking system and a simple thumb turn lock. The handles are adjustable for arm length and hand position and can be rotated inwards in 2 steps of 22.5 degrees each to aid posture alignment. The brakes have a simple “pull to brake” and “push to park” function. It features a textile seat, and a locking system that stops it reopening once unfolded when being loaded or moved.

The 11" & 8" smooth-running wheels are equipped with 2 brakes for optimum speed control that are also lockable for parking. The detachable bag has a load capacity of 5kgs and can also be used as a convenient shopping bag.

Specifications	MOT62
Overall Width	71cm
Overall Length	77cm
Product Weight	12.4kgs
Width Folded	31cm
Seat Height	62cm
Seat Width	56cm
Height of Arm Supports	87cm – 110cm
Max. User Weight	200kgs
Wheels	285mm/235mm x 44mm
Max. Bag Load	5kgs

OPTIONAL ACCESSORIES



WATERPROOF BAG



ONE HANDED BRAKE



SPEED REDUCTION WHEELS

Sold and maintained by: