




TIGER

FOREARM WALKER

TIGER FOREARM WALKER

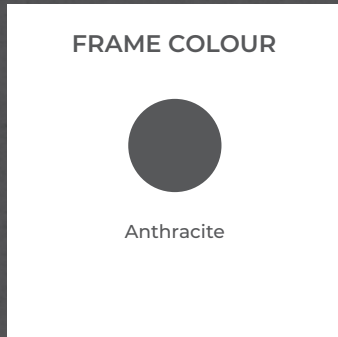

The Mobilex Tiger forearm walker is designed for use indoors where the ground surface is flat or smooth. It is perfect for users that require the additional support of the forearm rests, either as a result of weak wrist strength or rehabilitation purposes.

The forearm rests are height-adjustable, with an easy to read numerical marking system and a simple thumb turn lock. The handles are adjustable for arm length, hand and wrist position to ensure maximum comfort.

The large 200mm (8") smooth-running wheels have soft tyres for an enhanced ride. The rear wheels are equipped with brakes for optimum speed and movement control and the brake levers have a simple "pull to brake" and "push to park" function. As an option, a speed reduction system can be fitted for extra control.

Specifications	MT62
Overall Width	62cm
Overall Length	66cm
Product Weight	9.9kgs
Width Folded	22.5cm
Seat Height	62cm
Seat Width	46cm
Height of Arm Supports	95cm – 120cm
Max. User Weight	150kgs
Wheels	200mm x 35mm
Max. Bag Load	5kgs

OPTIONAL ACCESSORIES



NON-SLIP TRAY



ONE HANDED BRAKE



SPEED REDUCTION WHEELS

Sold and maintained by: