




TIGER PLUS

FOREARM WALKER

TIGER PLUS FOREARM WALKER

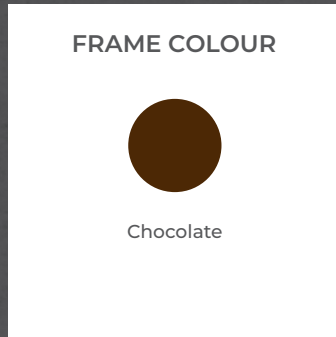

The Mobilex Tiger Plus forearm walker is designed for use indoors where the ground surface is flat or smooth. It is perfect for users that require the additional support of the forearm rests, either as a result of lack of wrist strength/support or rehabilitation purposes. Although very robust, it is still relatively lightweight at just 9.5kgs and offers a high degree of stability and ease of movement.

The forearm rests are height-adjustable, with an easy to read numerical marking system and a simple thumb turn lock. The handles are adjustable for arm length and hand position and can be rotated inwards in 2 steps of 22.5 degrees each to aid posture alignment. The brakes have a simple “pull to brake” and “push to park” function. It features a textile seat, and a locking system that stops it reopening once unfolded when being loaded or moved.

The 8" smooth-running wheels are equipped with 2 brakes for optimum speed control that are also lockable for parking. The detachable bag has a load capacity of 5kgs and can also be used as a convenient shopping bag.

Specifications	MT62P
Overall Width	62cm
Overall Length	66cm
Product Weight	9.9kgs
Width Folded	22.5cm
Seat Height	62cm
Seat Width	46cm
Height of Arm Supports	95cm – 120cm
Max. User Weight	150kgs
Wheels	200mm x 35mm
Max. Bag Load	5kgs

OPTIONAL ACCESSORIES



NON-SLIP TRAY



ONE HANDED BRAKE



SPEED REDUCTION WHEELS

Sold and maintained by: