



Confidence in posture solutions

# Azalea®



Offers excellent postural support  
and comfort with many adjustment  
possibilities!



The **Rea Azalea** is the natural choice for managing all aspects of postural support and comfort. It is available in three chassis sizes with options of both transit wheels for attendant propulsion, or self propel wheels for the moderately active individual. The seating system is versatile in its wide choice of support options and in the adjustments it offers. Making the **Rea Azalea** the market leader in tilt in space wheelchairs.

## Tilt & recline makes life comfortable

The **Rea Azalea** offers 25° of seat tilt and 30° of backrest recline to accommodate a range of positioning requirements. Both tilt and recline are easily adjusted by attendant operated gas pistons, supplied as standard. As an option, there is also an electric version to give the client independence over their own position and comfort.

Pivot points on the chair closely match the natural pivot points of the person, ensuring synchronised movements between the client and the chair, greatly reducing the risk of sheer during tilt and recline. For additional safety and confidence for the client, the unique Dual Stability System (DSS) of the **Azalea** incorporates weight shifting when tilting and reclining, to ensure stability is maintained.

## Features and Options



### Flex 3, backrest plate

Sturdy and supportive, with both width and height adjustments. Available with a range of cushions.



### Flex 3, tension adjustable backrest

To accommodate more complex spinal shapes, the tension adjustable backrest offers multiple adjustments, suitable for male and female body shapes of all sizes.



### Wide range of options

In addition to the multiple adjustments in the chair the Azalea boasts a large assortment of options and accessories, ensuring there is a solution for a wide range of postural challenges.

## Technical data


For more comprehensive pre-sales information about this product, including the product's user manual, please visit [www.invacare.co.uk](http://www.invacare.co.uk)

	Seat width	Seat depth	Seat height	Backrest height	Armrest height	Length of footrests	Seat angle	Backrest angle
<b>Rea Azalea Transit &amp; SP</b>	390-440/ 440-490/ 490-540 (390-490/440-540/ 490-590 mm)	430 - 500 mm	400/450 mm	Backrest plate: 560-790 mm Tension adjustable backrest: 540-715 mm	240 - 360 mm	330 - 500 mm	-1° to 23° (+1° to 25°)	0°-30°

	Driving unit width	Total height	Total length inc. Footrests	Total weight	Max. User weight	Transport weight	Crash tested	Maximum slope for use of parking brake
<b>Rea Azalea Transit &amp; SP</b>	SW + 220 mm (SW + 230 mm)	900-1250 mm	950-1020 mm	34 kg (SW 440 mm)	Max 135 kg	16,5 kg <sup>1</sup>	ISO 7176-19	7°

1. Transport weight increases with 4,2 Kg when the Power tilt & Recline option is used

## Frame colour

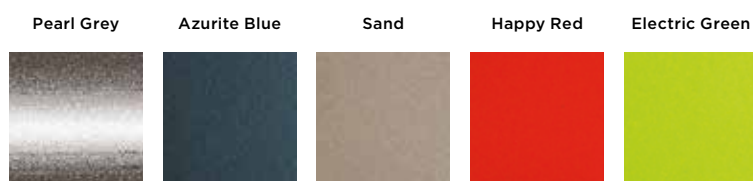

## Upholstery



Invacare Limited  
Pencoed Technology Park - Pencoed  
Bridgend CF35 5AQ - United Kingdom  
Tel: +44 1 656 776 222  
Fax: +44 1656 776 220  
E-mail: [uk@invacare.com](mailto:uk@invacare.com)  
Sales Order E-mail: [ordersuk@invacare.com](mailto:ordersuk@invacare.com)  
[www.invacare.co.uk](http://www.invacare.co.uk)

Invacare Ireland Ltd  
Unit 5 - Seatown Business Campus  
Seatown Road - SWORDS - County Dublin - Ireland  
Tel. +353 1 8107084  
Fax +353 1 8107085  
Email: [Ireland@invacare.com](mailto:Ireland@invacare.com)  
[www.invacare.ie](http://www.invacare.ie)



Yes, you can.®