

# J3 CARBON

## FEEL THE DIFFERENCE

# IT'S SO LIGHT...

---

Carbon fibre is lightweight and strong – the ideal material for an active backrest. Offering the same strength as aluminium (but 60% lighter), the J3 Carbon backrest reduces the weight you have to push. Day-by-day, you'll feel the difference!



## ... AND STYLISH

---

The J3 Carbon backrest will impress you with its modern design and clean finish, ensuring that the backrest is a perfect match for you and your wheelchair.



# THE SUPPORT YOU NEED TO BE ACTIVE

---

Designed to provide an ultra lightweight and stylish backrest for active users, the J3 Carbon backrest provides firm stabilisation of the pelvis and lower spine to encourage an optimal spinal curve. This stability is needed for upper body movement to ensure you can maintain an active lifestyle. Apart from being exceptionally comfortable, it can also help to reduce back pain.

The force that's generated when propelling your wheelchair is transferred more efficiently through a solid back into the frame. That means when pushing, you're generating greater movement with less effort.

**Just try it...**



# ADJUSTABLE TO YOUR NEEDS

---

## Choose your back...!

Different needs demand different solutions. The four back heights (from 17 - 38cm) focus on the active user with a need for low to moderate posterior trunk support. Mild lateral trunk support is provided by the 5cm contour depth.




→ **Active contour**

Mild lateral trunk support through 5cm contour depth.



→ **Different support levels**

Always the right support (four heights from 17 - 38cm).

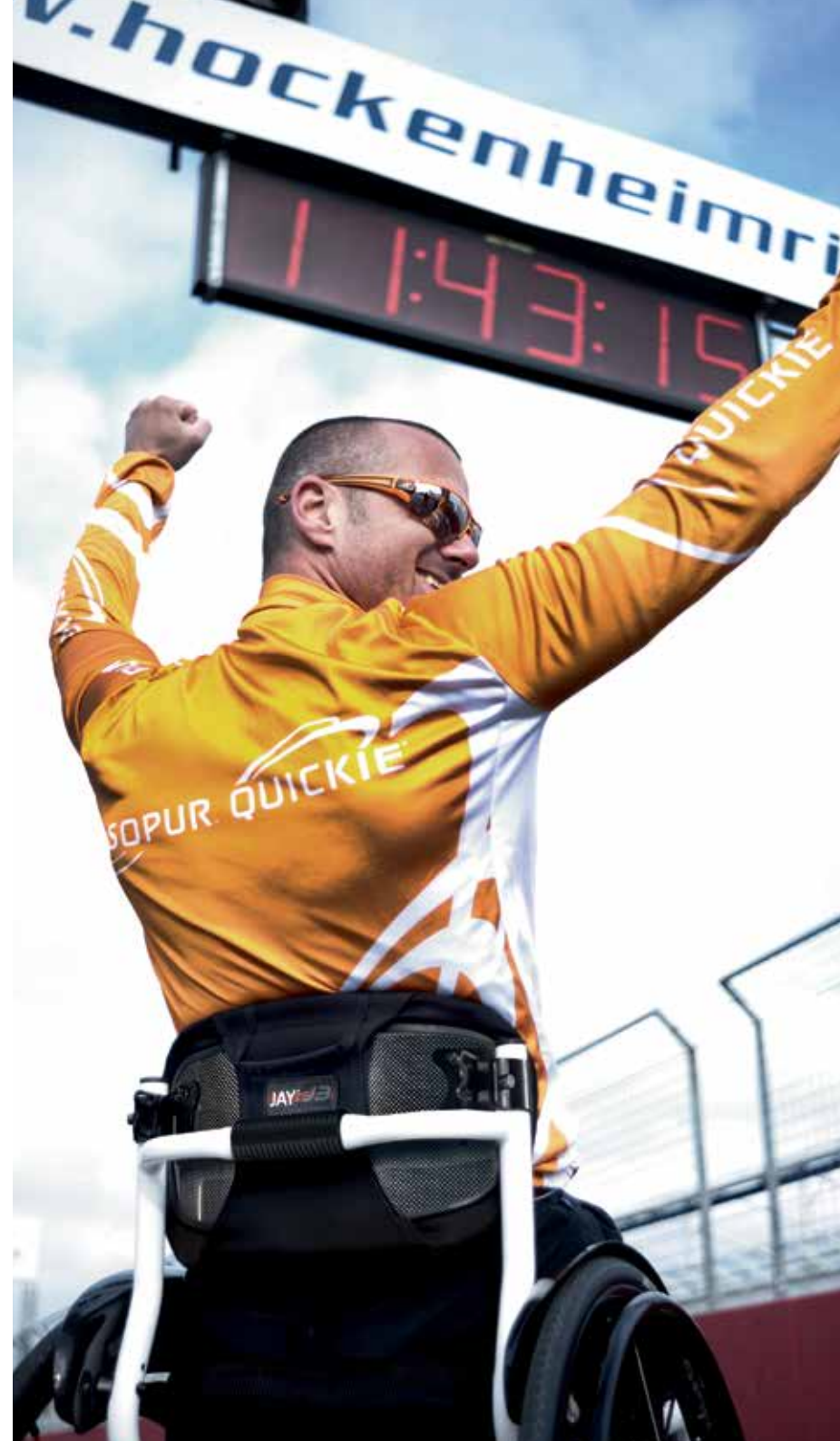

# CHOOSE YOUR MOUNTING SYSTEM

---

The J3 Carbon backrest is available with two different mounting systems. The fixed system offers rigid wheelchair users the lightest possible mounting solution. Alternatively, the quick-release system allows the backrest to be easily removed, ideal for folding wheelchairs. Both mounting systems offer the same angle, depth, height and width adjustments to adapt to your individual needs.



## The quick release hardware

Ideal for folding wheelchairs.

---



## The fixed mounting system

The minimum weight solution - ideal for rigid wheelchairs.

---



# LOOK COOL, STAY COOL


---

## Stay cool and comfortable with the microclimatic cover and soft foam underlay

The J3 Carbon backrest's cover uses 3DX™ spacer fabric, which vents both heat and moisture and increases air flow through the back. Machine washable at 60°C and quick to dry, the cover overlaps the shell to provide you with extra protection. You can even adapt the shape of the cover to your individual needs with lumbar support pads.

### → J3 Carbon cover

The 3DX™ technology offers a comfortable microclimate.

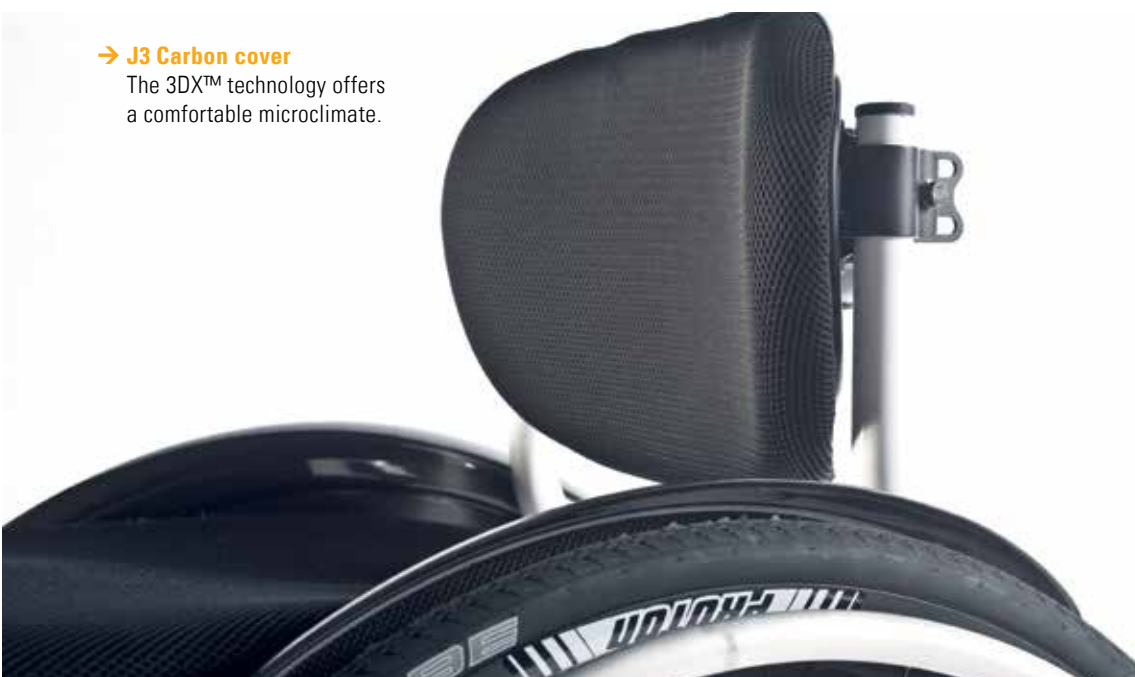

## Where are my keys?

You now have somewhere to store those easy to misplace items with the new accessory bag. Not only will it keep all of your valuables in one place, but it also smartly covers the gap between your back and cushion. You'll wonder what you ever did without it.

Accessory bag - useful and stylish.



# TECHNICAL DATA

<b>Maximum user weight:</b>	125kg
<b>Backrest width:</b>	32 / 36 / 41 / 46cm
<b>Backrest height:</b>	17 / 24 / 30 / 38cm
<b>Contour depth:</b>	5cm
<b>Width adjustment:</b>	5cm
<b>Depth adjustment:</b>	4cm
<b>Back angle:</b>	± 20°
<b>Backrest weight:</b>	From 690g (32cm wide backrest with quick release hardware)
<b>Hardware:</b>	Fixed and quick-release versions
<b>Accessories:</b>	Accessories bag and lumbar supports

Find out more today!

Visit [www.SunriseMedical.co.uk](http://www.SunriseMedical.co.uk) for further information on product options, or contact your local dealer or representative.



**Get SOCIAL with QUICKIE and JAY**

Keep up-to-date with the latest news on Quickie and JAY products through our Facebook page and YouTube channel:



[www.facebook.com/QuickieWheelchairsUK](http://www.facebook.com/QuickieWheelchairsUK)

[www.youtube.com/quickiewheelchairseu](http://www.youtube.com/quickiewheelchairseu)

## Sunrise Medical Limited

Thorns Road, Brierley Hill, West Midlands, DY5 2LD

Tel: +44 (0) 845 605 66 88 Fax: +44 (0) 845 605 66 89

[www.SunriseMedical.co.uk](http://www.SunriseMedical.co.uk)

