



JAY LITE

LIGHTWEIGHT.
BREATHABLE.
COMFORTABLE.



ULTRA LIGHTWEIGHT...

The perfect complement to an ultra lightweight wheelchair, such as the Quickie Helium and Quickie Xenon, the JAY Lite is the lightest cushion in its class.

Ultra lightweight from 660g, and offering increased skin protection, stability and comfort, the JAY Lite is the ideal cushion for the active user.




... AND BREATHABLE

Always active and on the go, you will generate heat underneath your buttocks and thighs that can lead to skin breakdown.

The JAY Lite uses innovative cushion technology to promote air-flow under the user whilst sat in their wheelchair.

The 3DX™ microclimatic cover and Olotex “airflow” base vents both heat and moisture through the cushion’s materials. This, added to the cushions’ styling, ensures users stay cool and look cool!



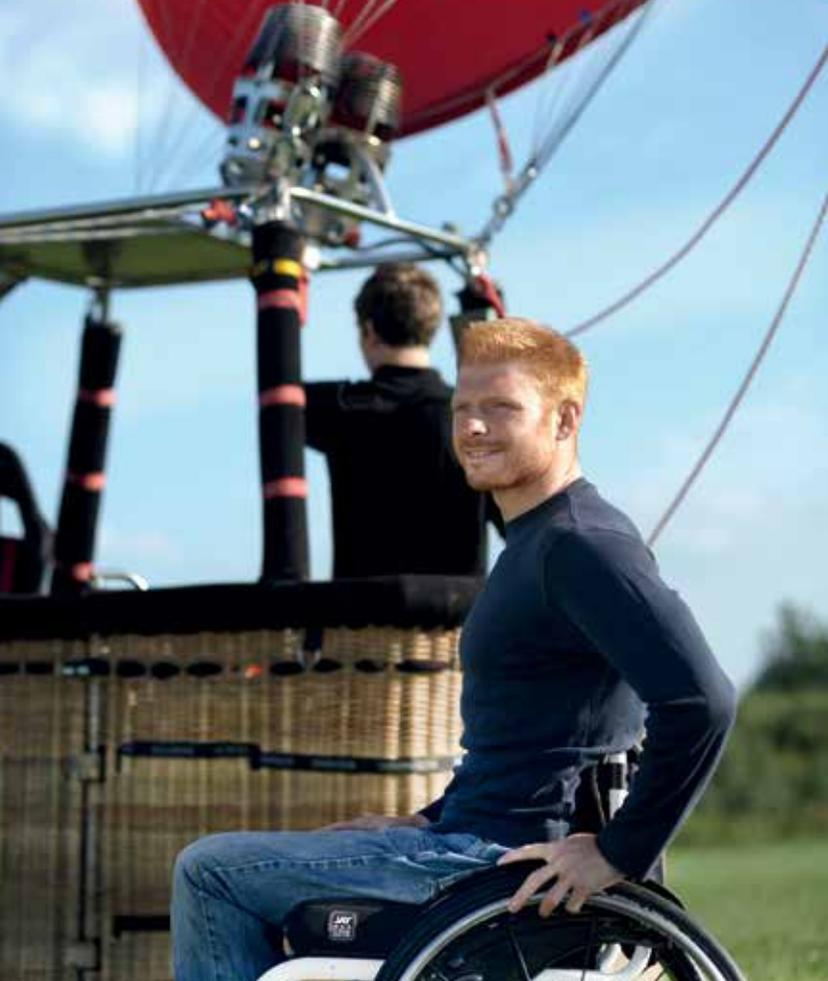

SUPERIOR POSITIONING

The JAY Lite cushion provides increased skin protection, stability and comfort. To minimise the risk of pressure sores, the anthropometrically designed Optiwell pelvic loading area (PLA) effectively reduces loading on the ischial tuberosities without overloading the hip bones.

Apart from increasing breathability, the lightweight reticulated foam insert features a 'scoop-out' which helps to redistribute pressure under the seat bones. Plus, the pre-contoured foam base provides the support you need to keep active and comfortable!

Pelvic loading area


Redistributes weight to reduce the risk of pressure sores

Sacral notch

Protects the coccyx area

Ischial well cutouts

Reduces pressure under the ischial tuberosities to protect skin integrity

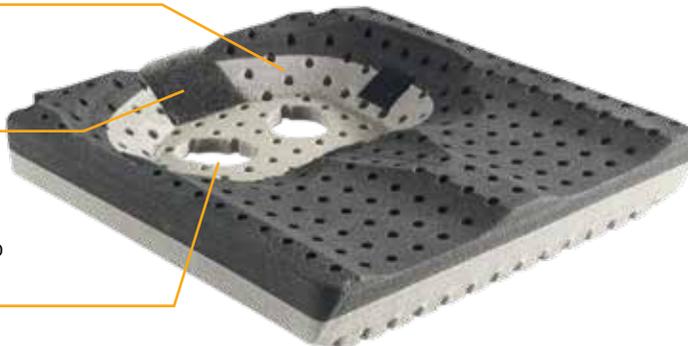


Reticulated foam insert

Superior breathability, washable and removable

Pre-contoured base

Increased stability and leg support

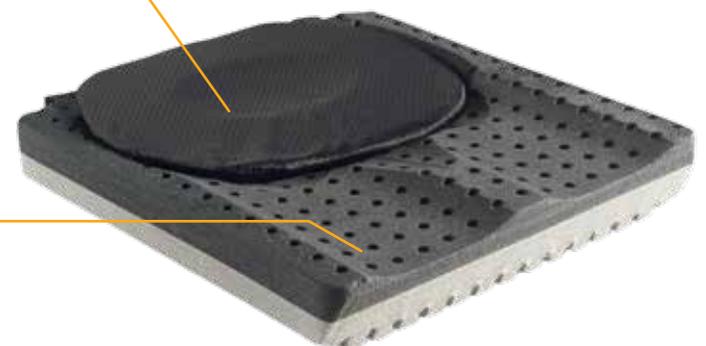



ONE CUSHION WITH THE CHOICE OF TWO DIFFERENT BASE CONTOURS.

The JAY Lite cushion offers two different contoured base versions. The standard base is ideal for most users who require comfort and pressure redistribution. The standard low-profile contour allows for easier transfers in and out of the wheelchair. Alternatively, for those who require a higher level of postural stability, the JAY Lite P version provides deeper leg troughs which increases the overall stability of the thighs and pelvis.

→ **JAY LITE – low profile contour**
Provides mild leg support and allows easy transfers.



→ **JAY LITE P – positioning contour**
More leg positioning through optional lateral and medial thigh support for increased overall stability.



INNOVATIVE AND FUNCTIONAL DESIGN

The JAY Lite's design meets the lifestyle needs of an active user - lightweight, low profile, breathable, pressure relieving and stable. The bevelled front edge prevents rubbing against the back of the legs, eliminating the need for using a shorter length. The corner cut at the lower sides of the cushion also reduces the load on the trochanter (and reduces the seat-to-floor height). The custom carrying handle offers an easy grip and the notches for back canes locks the cushion into place.



Grey rubber non-stick bottom
Breathable and anti-slip



Corner cut
To reduce load on the trochanters



Bevelled front edge
Prevents rubbing against the back of the legs

Custom carrying handle
For easy grip

Select your cushion insert.

The JAY Lite cushion offers two different well inserts.

The reticulated foam insert uses an air-flow cover and optimises breathable airflow. The entire insert is also machine washable at 60°C and quick drying.

The optional PU foam insert comes with an incontinence cover that's water-resistant and easy to wipe-off and clean.



→ Reticulated foam insert

With air flow cover.
Superior breathability,
machine washable. Central
scoop-out reduces pressure
under seat bones or ischials.



→ PU foam insert (optional)

With incontinence cover.
Insert cover is water-
resistant and easy to
wipe off and clean.

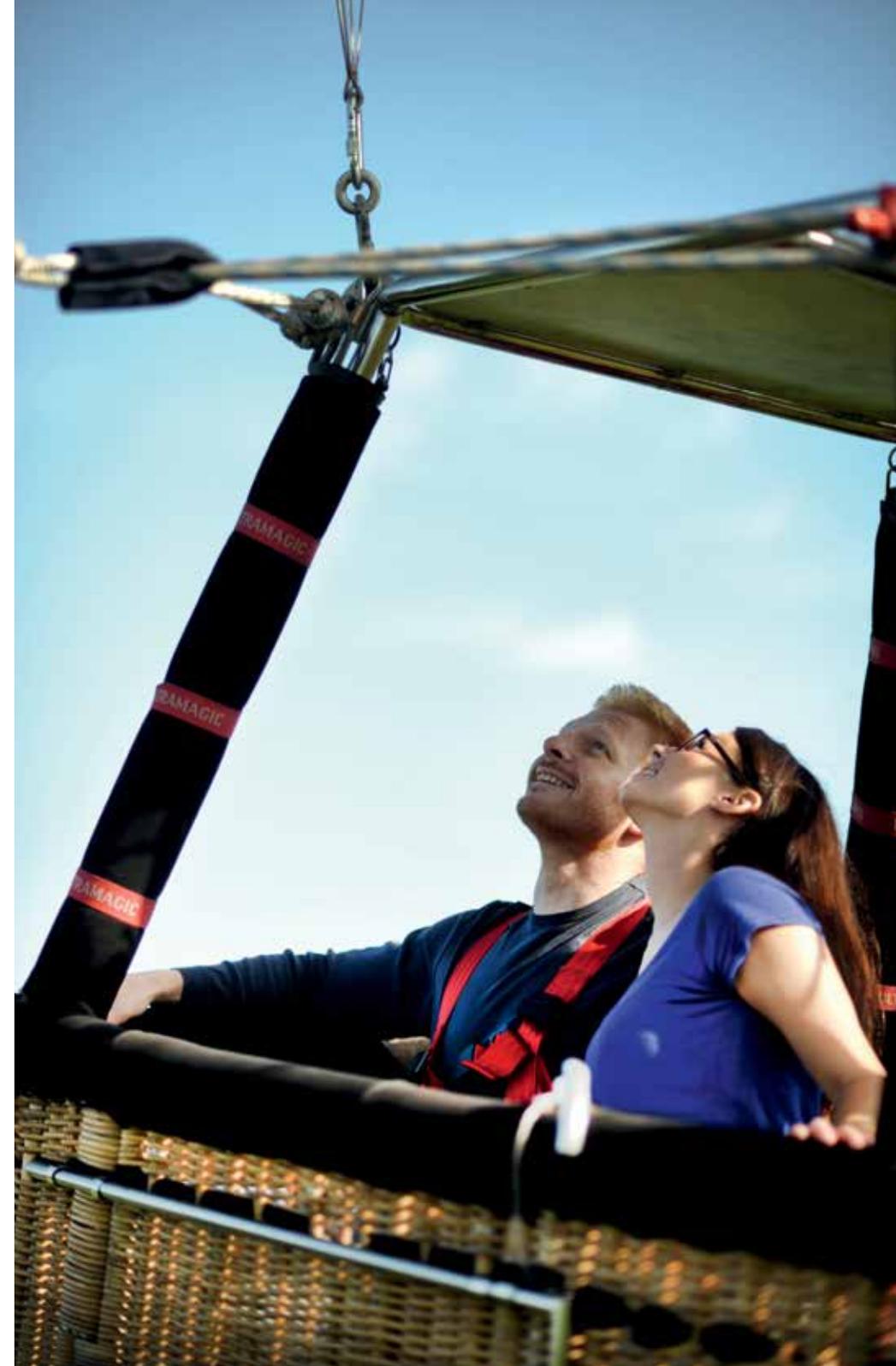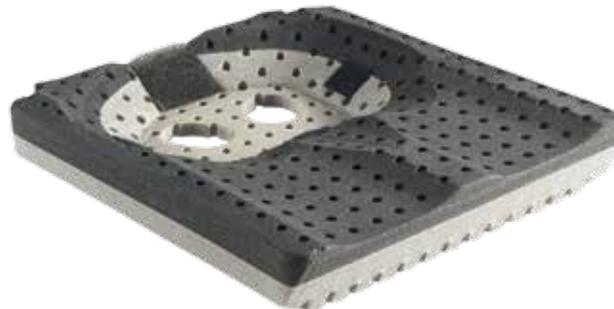

TECHNICAL DATA

Cushion width:	35cm – 50cm
Cushion depth:	35cm – 50cm
Cushion height:	8cm
User weight capacity:	150kg
Cushion weight:	660g - 1090g
Insert:	Reticulated foam insert with air flow cover (standard) & PU foam insert with incontinence cover (optional)
Outer cover:	Microclimatic cover
Base versions:	Standard version / Positioning version - contoured base for better leg positioning

Find out more today! Visit www.SunriseMedical.co.uk for further information on product options, or contact your local dealer or representative. For further information on the full specification, options and accessories, please refer to the order form. All information is subject to change without notification. Please consult Sunrise Medical with any queries you may have.

Get SOCIAL with QUICKIE and JAY

Keep up-to-date with the latest news on Quickie and JAY products through our Facebook page and YouTube channel:

 www.facebook.com/QuickieWheelchairsUK

 www.youtube.com/quickiewheelchairseu



PB_JAY_LITE_CUSHION_2014_REV1

Sunrise Medical Limited

Thorns Road, Brierley Hill, West Midlands, DY5 2LD
Tel: +44 (0) 845 605 66 88 Fax: +44 (0) 845 605 66 89

www.SunriseMedical.co.uk