



Life Family

Live life ...
... without limits



QUICKIE[®]
Live without limits.

Life outside the box

A complete fleet solution without compromising on performance...

- **Performance** - compact, lightweight with on-chair adjustments for fine tuning, Life delivers excellent rolling performance for active users
- **Lightweight** - the lightweight frame and new range of lightweight options enable easy transfers and lifting for active users
- **Unparalleled adjustability** - the widest range of easy on-chair adjustments to meet changing clinical needs, improving wheelchair skills or reissue
- **Stylish design** - a modern take on the classic frame design complemented by modern options borrowed from our high end products, Quickie Helium and Xenon

Life Rigid

for performance



Easy Prescription



Easy Adjustments

Life Folding

for transportation



Efficient Fleet Management: The fleet approach to active clients

The Quickie Life family provides a flexible and high-performance fleet solution, which can be perfectly configured to meet your active client's needs...

- ✓ Completely interchangeable parts between Life Folding and Rigid
- ✓ Modular frame allows growth and crossover of options to meet multiple client needs and recycling
- ✓ Widest range of adjustability built in to every chair
- ✓ Cross compatible spares and options
- ✓ **The result... reduced complexity and greater flexibility**



Easy Reissue



Easy Life

Performance...

Lightweight...

Energy Efficient...

Configurable...

Robust...



Excellent rolling performance via compact adjustable wheel base



Three stabiliser bars deliver rigidity for performance and longevity



Fixed front frame for ultimate rigidity and performance.



Swing-away option/Xenon lightweight hangers with easy to use trigger mechanism.



Small folding with a double locking system for easy lifting



Rear frame protection to protect the frame for transfers

... and Durability

Crash tested to ISO 7176-19 and passed:

- ✓ Rolling road **200,000 cycles**
- ✓ Drop test **6,666 cycles**
- ✓ Static stability - tested on a tilting table under maximum user weight
- ✓ Static load - including:
 - ✓ footrest downward load - 125kg required **155kg achieved (+20%)**
 - ✓ footrest upward load - 102kg required **108kg achieved (+6%)**
 - ✓ armrest downward load - 97kg required **107kg achieved (+11%)**
 - ✓ footrest impact 10kg pendulum from 51°
 - ✓ castor resistance 10kg pendulum from 51°
 - ✓ backrest impact 25kg pendulum from 30°

...the result... like our Quickie GPV, a high performance active chair, that will last the test of time



Unparalleled adjustability... making life easy



Adjustable Axle Plate

Built in adjustment for fine tuning performance or recycling for multiple clients

- 8cm rear seat height adjustment built into every chair
- Seat height adjustable in 7 positions
- 6cm COG adjustment (2.3 to 8.3cm)
- 1° to 4° camber adjustment requiring no new parts

Backrest Angle Adjustment

- -12° to +8° in increments of 4° - Life R
- -12° to +12° in increments of 4° - Life F
- Lightweight and simple to adjust
- Pivot point 6cm above the seat to match the hip angle
- Double locking system for easy transportation



Life R



Life F



Backrest Height Adjustment

- Tension adjustable backrest as standard (Helium upholstery)
- 25 - 45cm of back height adjustment across 4 options
- 5cm on-chair adjustment





Seat Upholstery Adjustment

- 36 - 50cm of seat depth adjustment via the canvas (replacement front frame required)



Castor Adjustment

- Easy castor fork angle adjustment to accommodate changes in seat to floor height
- Castor height adjustment via stem bolt and 4 positions on fork
- Set forward castor fork position to increase stability

Hangers

- Choice of **swing-away** legrests or **fixed front frame** to accommodate user preference
- Lower leg length adjustment available from 27 - 47cm dependent on options



Manual elevating legrest from Xfamily²



70° or 80° hanger, Xenon (shown) or HeliX² style



75° or 85° fixed front frame. 2cm inset option

Lightweight Options

Lightweight product for active users ...

... designed to enrich your life ... for life

The Life and Life R manual chairs offer a wide range of lightweight options enabling the chair to be customised in a way to suit everyone.

- » three arm cross brace (folding Life) and three towel bars (rigid Life R)
- » lightweight wheels and handrims
- » lightweight platform footrest - fixed front
- » lightweight aluminium hangers - swing-away (Xenon style)
- » NEW lightweight armrest and side guards (Helium style)
- » NEW lightweight anti-tips - **70% lighter!**

The modern frame design and lightweight options minimise energy and forces needed by the user for everyday activities, reducing the loads on shoulders, arms and hands.



NEW Tool height adjustable armrest
Lightweight modern looks with one pad for easy transfer and low transport weight



NEW Aluminium lightweight hangers
in 70° and 80°. 13% lighter than standard Quickie hangers



Lightweight anti-tip for optimal driving performance with 480g weight reduction.



Desk armrest with improved stability, integrated handle and improved armpad position (4cm forward)

Armrests and Side Guards



Tool height adjustable armrest

- Lightweight - only 730g each
- One pad for short and long position
- Easy transfer and low transport weight



Single post height adjustable armrest

- Robust with increased comfort
- Short or long pad choices
- Easy transfer and low transport weight



Composite side guard

- Detachable with low transport weight
- Ideal for side transfers
- Good clothing protection



Aluminium side guard

- Low weight
- Good clothing protection



Aluminium side guard with fender

- Higher safety with fender
- Good clothing protection
- Provides support for transfers

Footrests



Lightweight platform footrest

- Angle and depth adjustable
- Increased rigidity
- Foot positioning plates (Optional)



Platform footrest

- Angle and depth adjustable
- Easy transfers



Lightweight platform footplate

- Angle and depth adjustable
- Easy transfers (flip-up)
- Low weight



Divided footplates (Aluminium)

- Angle and depth adjustable
- Easy handling and transfer (flip-up)



Divided footplates (Composite)

- Angle and depth adjustable
- Easy handling and transfer (flip-up)
- Less weight impact

Colour options

The Quickie Life and Life R are available in the following fantastic colours, helping you to customise the chair.



White (01)
(RAL 9010)



Brilliant Silver (02)



Brilliant Lime (07)




Matt Marine (16)



Matt Red (23)



Black (29)
(RAL 9005)

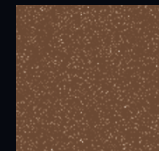

Dark Chocolate
Metallic (34)



Cyan Metallic (36)



Mars Red (37)



Matt Anthracite (42)

More adjustment per chair >>> Interchangeable parts including front frame >>> Easy growth >>> Easy re-issue

More Efficient Fleet Management

The Life Folding & Life Rigid use the same modular frame design allowing you to change seat depth from swing-away to fixed front or alter the frame type by simply changing the front of the frame.


Efficient Fleet Management: The fleet approach to active clients

With compatible parts between Life Folding and Life Rigid and the modular frame design means not only can you change options and accessories, but you can also change seat depth or front frame type (swing-away or fixed) by simply and cost effectively changing the front frame, as shown in the example to the right.

- ✓ Completely interchangeable parts between Life Folding and Rigid
- ✓ Modular frame allows growth and crossover of options to meet multiple client needs and recycling
- ✓ Widest range of adjustability built in to every chair
- ✓ Cross compatible spares and options


Ultimate growth with interchangeable frames

Step 1

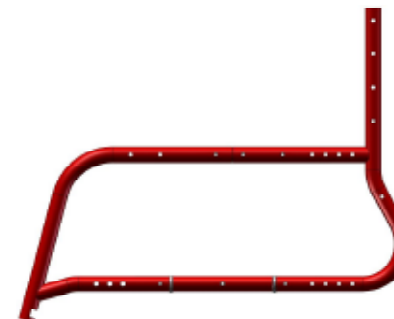


Example shows a fixed front frame

Step 2

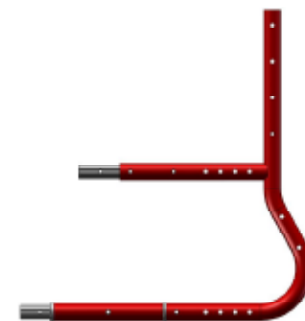

Remove front end from frame

Step 3

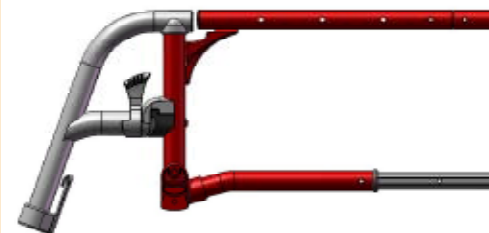

Connect new swing-away front end to rear frame

The complete approach to fleet management ... 4 chairs in 1... with great cross compatibility

Component / option	Folding		Rigid	
	Life SA	Life FF	Life R SA	Life R FF
Seat - Life	✓	✓		
Seat - Life R			✓	✓
Back sling	✓	✓	✓	✓
Back (standard)	✓	✓	✓	✓
Push handles	✓	✓	✓	✓
Composite hanger 70°	✓		✓	
Composite hanger 80°	✓		✓	
Aluminium hanger 70°	✓		✓	
Aluminium hanger 80°	✓		✓	
Elevating legrest	✓		✓	
Footplates	✓	✓	✓	✓
Side guards	✓	✓	✓	✓
Tool height adjustable armrest	✓	✓	✓	✓
Height adjustable armrest	✓	✓	✓	✓
Desk armrest	✓	✓		
Wheel locks, extension lever, castors, forks, rear wheels, axle plate, spoke guards, anti-tips, tip assist	✓	✓	✓	✓



Also shares common parts with the Quickie Helium, Xenon, Argon, Neon and HeliX² manual chairs extending the fleet management solution over even more products!



	Life fixed front	Life swing-away	Life R fixed front	Life R swing-away
Overall width	Seat width +20cm (8")	Seat width + 20cm (8")	Seat width + 20cm (8")	Seat width + 20cm (8")
Seat width	36 - 50cm (14 - 20")	36 - 50cm (14 - 20")	34 - 48cm (13 - 19")	34 - 48cm (13 - 19")
Seat width increments	2cm (0.7")	2cm (0.7")	2cm (0.7")	2cm (0.7")
Back height	25 - 47.5cm (10 - 19")	25 - 47.5cm (10 - 19")	25 - 45cm (10 - 18")	25 - 45cm (10 - 18")
Seat depth	36 - 50cm (14 - 20")	36 - 50cm (14 - 20")	36 - 50cm (14 - 20")	36 - 50cm (14 - 20")
Seat depth increments	2cm (0.7")	2cm (0.7")	2cm (0.7")	2cm (0.7")
Rear seat height range	40 - 50cm (16 - 20")	40 - 50cm (16 - 20")	36 - 45cm (14 - 18")	36 - 45cm (14 - 18")
Front seat height range	40 - 54cm (16 - 21")	40 - 54cm (16 - 21")	39 - 50cm (15 - 20")	39 - 50cm (15 - 20")
Back angle range	-12 to +12°	-12 to +12°	-12 to +8°	-12 to +8°
Frame front end angle	75°, 85°, 85° with 2cm inset	80° and 70°	75°, 85°, 85° with 2cm inset	80° and 70°
Individual frame colours	10	10	10	10
Individual sling colours	1	1	1	1
Fork and axle plate anodized colours	3	3	3	3
Total chair weight, from	11.5kg	12kg	10.7kg	11.3kg
Maximum user weight	125kg (19.6 stone)	125kg (19.6 stone)	125kg (19.6 stone)	125kg (19.6 stone)

Sunrise Medical
 Thorns Road
 Brierley Hill
 West Midlands
 DY5 2LD
 Tel: +44 (0) 845 605 66 88
 Fax: +44 (0) 845 605 66 89
www.SunriseMedical.co.uk

www.SunriseMedical.co.uk

