

ottobock.

Start Multi.

The new standard in care.

The *Start Multi* is an innovative multifunctional wheelchair that offers maximum comfort and adaptability. It is available in two versions (12" and 24" rear wheels) and three sizes (small, medium and large).

The seat frame allows tilting from -5° to 30° , while the back support can be adjusted continuously from 92° to 128° .

Ergonomic seat cushion and back support set a new standard for comfort. The elevating leg supports flexibly cater to individual needs. Additional features are adjustable head support, adjustable arm support, a pelvic belt for safety and proper positioning, and two anti-tippers for maximum safety during use. It is easy to operate, both for users and caregivers.

Our multifunctional wheelchair has been designed to meet the challenges of caregiving while fulfilling the highest standards of comfort and adaptability.



Start Multi.

The new standard in care.



Seat Angle Tilt

Tilting via gas spring from -5° to 30° is a standard feature for the Start Multi. Release with hand lever.




Adjustable Back Recline

Back recline continuously adjustable via gas spring from 92° to 128°. Easy release with hand lever.



Elevating Leg Supports

Five predefined positions of elevating leg support for symmetrical positioning. The new leg support is easy to release & adjust.

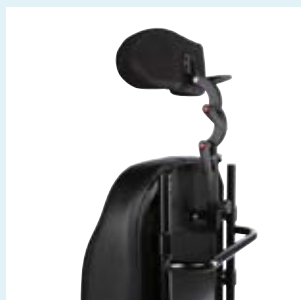

Head support

Head pads with high-resilience foam and breathable cover. Articulated mounting kit for height and depth adjustments.



Height adjustable arm support

The arm support is height-adjustable, easily adjusted by pressing the side button, and the injected foam cushion can be adjusted to different depths.

Other features:

- Abductor
- Anti-tipper
- Adjustable foot plate
- Telescopic push handles
- Wipeable seat and back cushion cover
- Adjustable center of gravity (different positions for the rear wheel)
- Adjustable caster wheel

Technical data

Artikelnummer	Size	Seat width (cm)	Rear Wheel
482M166=00001*	Small	39 (adj. to 45cm)	12"
482M166=00002*	Medium	44 (adj. to 50cm)	12"
482M166=00003*	Large	49 (adj. to 55cm)	12"
482M166=00004*	Small	39 (adj. to 45cm)	24"
482M166=00005*	Medium	44 (adj. to 50cm)	24"
482M166=00006*	Large	49 (adj. to 55cm)	24"