

Permobil

M Corpus[®] VS



permobil

Designed with the user in mind

The M Corpus® VS power wheelchair brings a new dimension to power standing from Permobil, with the maneuverability of a mid-wheel drive base.

Stand without compromise

The M Corpus VS combines a compact base and suspension that provides stability, traction, and ride comfort for a best-in-class power standing experience with a mid-wheel base. Designed with the user in mind, the M Corpus VS allows easy adjustment of the standing sequence & positioning to your individual needs in real time.

From face-to-face interactions to improved health outcomes, the M Corpus VS empowers you to stand with confidence, while our light, swing-away knee supports offer improved comfort and flexible on-chair storage for greater independence.

Designed for independence

With a lightweight design, our swing-away knee supports can be moved and handled independently using the easy to reach release mechanism even with limited hand function. With two pad sizes, choose the option that best matches you and easily fine tune the pad placement as needed with the integrated BOA® dial* for optimal fit and comfort to your needs.

Programming simplified


QuickConfig is a new innovative app designed for service technicians, clinicians, and providers that offers a simpler way to configure parameters and an intuitive interface to set up M Corpus VS with your preferred driving and seating characteristics and standing sequence.



M Corpus® VS Specifications

Max. speed	6 and 10 km/h
Range (G24: 85 Ah batteries)	40 km*
Max. user weight	120 kg
Obstacle climbing (with run up)	50 mm (60 mm)
Min. turning diameter	1310 mm
Turning in corridor	1000 mm
Suspension	Duallink
Max. standing angle	70°
R-net electronics Power Platform	120 A
Width incl. seat	650–790 mm
Height	1090–1170 mm
Min. transport length	880 mm
Min. transport height (folded backrest)	825 mm
Weight incl. batteries 85 Ah	201 kg
Weight of batteries 85 Ah	2x23 kg
Seat height (with electric seat elevator)	450–750 mm
Seat depth	370–570 mm (by 25 mm)
Seat width	420–570 mm (by 50 mm)
Backrest width	360/410/460/510 mm
Backrest height	470, 545-670 mm (by 25 mm)
Armrest height	215-330 mm
Distance between armrests	380–480/480–580 mm

*Driving conditions and choice of battery type. Measurements are approximates and subject to change without notice.

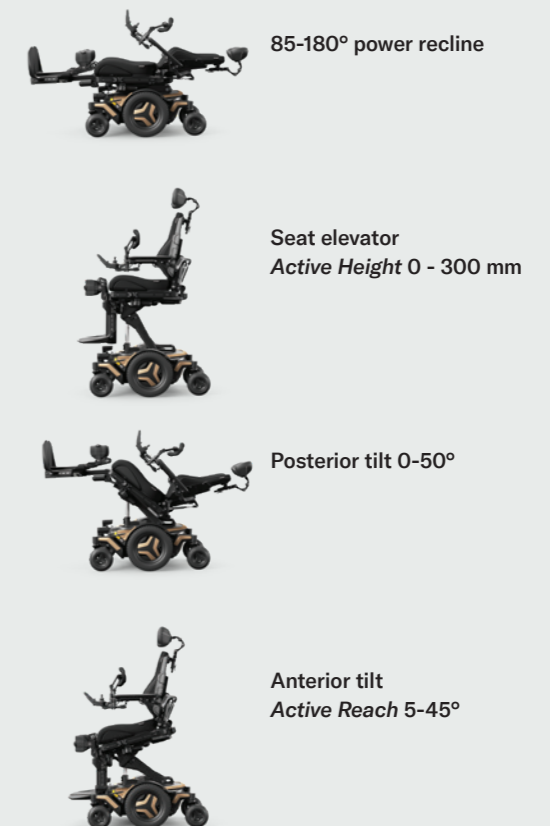

*BOA® Fit System delivers micro-adjustable, precision fit in a simple to use dial



Standing supports
face-to face interactions
and enables improved
psycho-social outcomes.

Permobil

M Corpus[®] VS

With the M Corpus VS, we've combined the benefits of our individualised standing functionality with mid-wheel drive base maneuverability. With a shorter overall length and smaller turning radius, experience improved accessibility to small and confined spaces.

permobil

customer_support_uk@permobil.com | permobil.com



GR3965-0.2UK230601