

Q300E

The compact powerchair that
doesn't downsize on adventure.



Q300 R



Q300 R Mini

DON'T DOWNSIZE ON ADVENTURE BECAUSE YOU WANT A COMPACT POWERCHAIR.

That's why with the **Q300 R**, you're not expected to. Indoors, it's highly manoeuvrable - super narrow and nimble in the most-crowded of situations. Yet, venture outside and it boasts the fierce traction, motor technology and drive performance to rival a much larger powerchair.

And it's our rear-wheel design that makes all the difference. A meticulously tested and refined base that provides the perfect balance between stability and size. It's also a base that can position your weight to maximise grip and conquer those big wide-open spaces. Prefer indoor agility over performance? Choose the **Q300 R Mini** and enjoy a slimmer **550 mm** base and 12" drive wheels for the tightest turning and manoeuvring.

Whichever you choose, its sophisticated **SEDEO LITE** seating system provides the same clinical support for the most challenging of conditions. Available with a host of power seat options and striking shroud and accent colours to dazzle, the **Q300 R** looks (and feels) amazing inside and out.



Scan the QR code with your
phone's camera to watch video.



ALL-ROUND PERFORMANCE

THOUGHT ALL REAR-WHEEL CHAIRS DROVE THE SAME? THINK AGAIN.

Take the Q300 R for a spin and you'll notice the difference. The way it grips the ground. The superb handling. But how does it achieve what others can't?

It's all thanks to the Q300 R's optimised base. By carefully testing and refining hundreds of precision adjustments to the width, length, wheel, and seat position, we found the smallest base geometry without compromising stability or traction.

Because the Q300 R is so stable, it can also support a wider range of centre-of-gravity adjustment. So, your seat position can be moved further forward or back to better balance your weight and increase traction and drive performance. Plus, you also get a compact chair that performs just as good indoors too!



Q300 R



Q300 R MINI

OUTDOOR PERFORMANCE

GET TO GRIPS WITH GREATER TRACTION.

With a precision-measured base that's geared for outdoor performance, the Q300 R has all the serious grip and ride performance you need to handle a busy city centre or off the beaten track.



KERB CLIMBING

No drop kerb in sight? Thankfully, the Q300 R with its boot-style kerb climber can smoothly handle kerbs of up to **100 mm**. Now you can say goodbye to taking the long way around.



14" DRIVE WHEELS

Looking for pure outdoor performance? 14" Drive wheels provide greater traction and stability over the harshest of terrain. They'll also keep your Q300 R slim at **580 mm**.



32 km

RANGE

Powered by (optional) punchy **63 Ah batteries**, you'll get everywhere you want to go without having to rush back with a range of up to **32km (20 miles)**.



REAR SUSPENSION

Providing responsive cushioning when driving over bumpy or uneven terrain, the Q300 R's rear suspension reduces the amount of shock and vibration that passes through your powerchair. So, you can sit-back and enjoy a soft and smooth ride.

KEEP ON TRACK WITH GYRO CONTROL

As there's fewer castor wheels to spin back into direction when turning, the Q300 R naturally drives in a straight line at higher speeds and covers ground more quickly without needing to oversteer. But when you do encounter a non-level surface, the optional Gyro Control will automatically make frequent steering corrections to keep you moving in a straight line - beneficial if you're using speciality controls or have poor hand co-ordination.



INDOOR MANOEUVRABILITY

AT HOME, INSIDE AND OUT.

From tight corridors, busy offices to tiny homes, the Q300 R's slim 580 mm base takes everything in its stride. With superb handling and intuitive rear-wheel drive, you'll find it incredibly easy to manoeuvre.

430 mm

SEAT-TO-FLOOR HEIGHT

Fit comfortably under tables or when driving your powerchair into your wheelchair adapted vehicle with the Q300 R's ultra-low seat-to-floor height that starts from **430 mm (17")**.



1800 mm

TURNING CIRCLE

1100 mm

BASE LENGTH

Ultra-short in length, the Q300 R fits in-and-out of the smallest spaces, so you no longer have to worry about taking up too much space behind.



GET EVEN MORE COMPACT WITH

Q300 R Mini

Looking for tighter manoeuvrability? Packing 12" drive wheels and centre mount as standard, the Q300 R Mini's ultra-slim 550 mm base can go where others can't. Plus, it's light enough (from 124 kg with 50 Ah batteries) to be compatible with most floor lifts too.



1060 mm
BASE LENGTH



550 mm
BASE WIDTH



GET CLOSER WITH THE CENTRE MOUNT LEGREST


With the centre-mount legrest, not only can you get closer to objects (like tables or serving counters), but you'll get the tightest turns too!

1700 mm
TURNING CIRCLE

POWER SEATING

SEDEO LITE, THE BEST SEAT IN THE HOUSE.

Combining an adjustable and configurable seat with the broadest choice of support and positioning options, the SEDEO LITE is the medium-to-high clinical seating solution that's designed to be tailored to your changing needs (and all with just one tool).

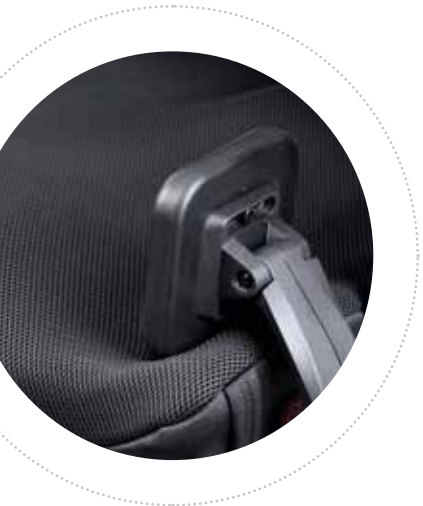

INTEGRATED SEAT RAIL

The integrated seat rail offers many mounting options for positioning elements like lateral thoracic, thigh and knee supports. All you need to do is attach and slide into position.



UNIVERSAL HEADREST

Compatible with **SEDEO** and **WHITMYER** headrests, this mount will help to address a range of head-support needs (from basic to multiple-point control).



LENGTH COMPENSATING, SWING-AWAY LEGRESTS

Adjustable in length to achieve the best fit, these legrests swing away for easier seat transfers. Also available as a power elevating and centre mount versions.

SEDEO PRO FLEXIBLE BACK

Upgrade to this unique backrest that effortlessly combines adjustable full rear back and side trunk support with the highest level of immersion that must be felt to be believed!

SHAPES TO YOUR SPINE

Unlike solid backs, the Flexible Back wraps around and contours to your spine's shape to provide complete contact, comfort, and stability. Because it wraps around you, the integrated laterals hug your sides, causing the foam to mould to your anatomical curves and keeping your trunk stable and aligned.



It's so immersive, even moderate-to-severe kyphotic postures or spinal scoliosis can be accommodated.



DO YOU HAVE SPECIFIC SEATING REQUIREMENTS?

The SEDEO LITE is also compatible with **J3 Back**. Suitable for moderate-to-high positioning, the J3 Back offers the widest range of sizes, shapes, and options to precisely fit and tailor your seat just to you.

POWER OPTIONS

POWER-UP YOUR POWERCHAIR

The Q300 R doesn't just have **BIG** chair performance, it's also loaded with **BIG** features to meet your clinical and lifestyle needs.



POWER ELEVATE AND TILT


Reach new heights with **300 mm** of seat lift and get to everyday objects that would otherwise be out of grasp. Transferring is easier too thanks to being in a higher position to stand from. Plus, with **30°** of tilt, you can sit back and enjoy greater relaxation and pressure distribution.





R-NET ADVANCED JOYSTICK

Upgrade your on-board technology with this large, easy-to-read HD display and simple to use paddle switches. Its built-in Bluetooth and i-device technology also allows you to connect to your smart devices all from the comfort of your Q300 R.

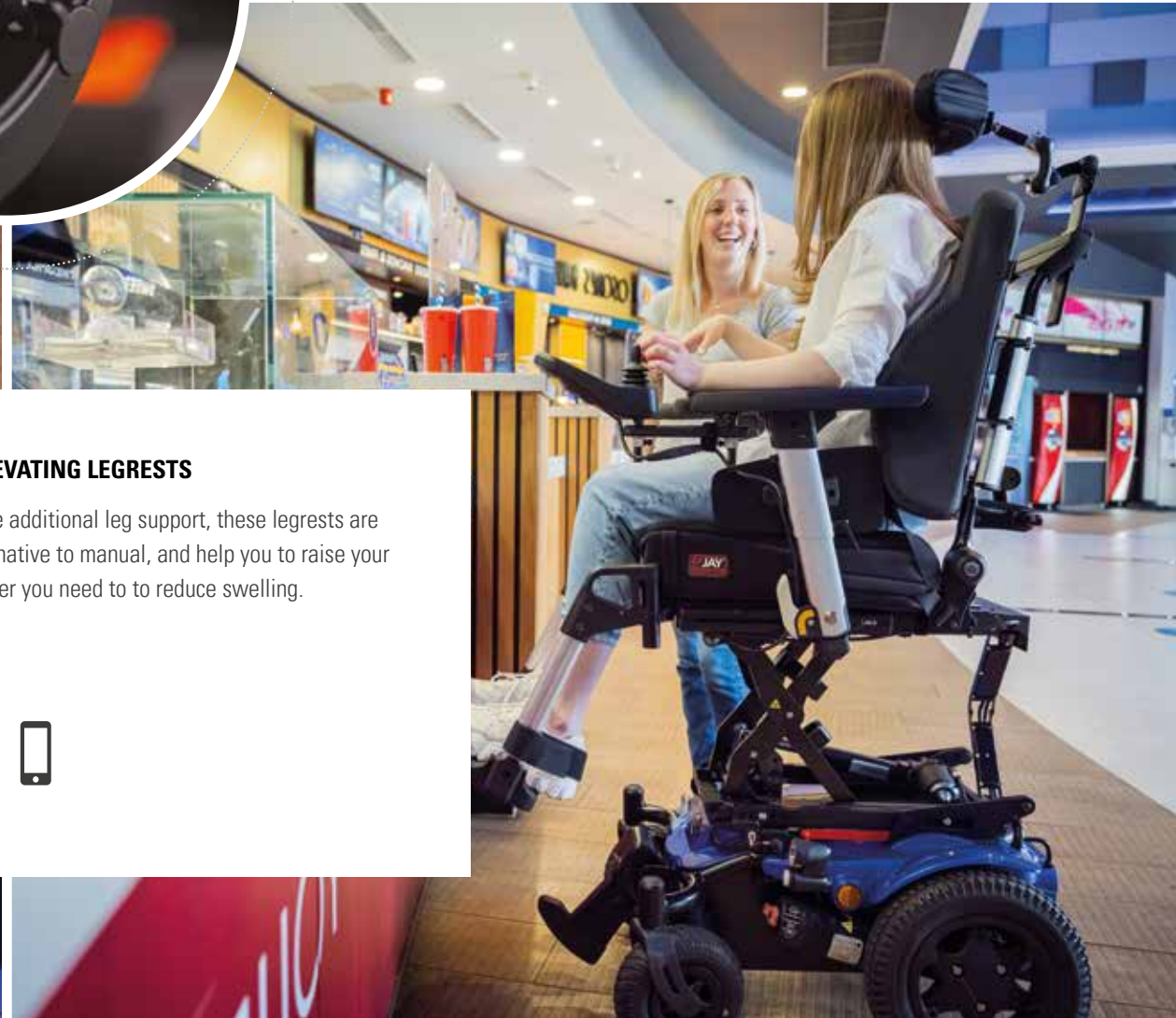

120° POWER RECLINE

Need to relieve pressure and re-distribute your weight throughout the day? Power recline enables you to adjust your seat-to-back angle in a similar manner to how a seat would recline in a car.



POWER ELEVATING LEGRESTS

If you require additional leg support, these legrests are a great alternative to manual, and help you to raise your legs whenever you need to to reduce swelling.





C-ME

ELEVATE AND DRIVE WITH C-ME

Unfortunately, life isn't designed around the height of your powerchair. That's why C-Me allows you to elevate up to 300 mm and drive up to 4.8 kph all at the same time! Now you can keep up with your friends at eye-level and at walking speed. Reach those heights that were out of grasp. And always be seen instead of overlooked.



SUSPENSION EVEN WHEN ELEVATED

With C-Me, the suspension on your powerchair remains active even when you're elevated. So, you get the same smooth and stable ride across a variety of terrain regardless of seat level.



ELEVATE WHILST DRIVING WITH ASSIGNABLE BUTTONS

Fed up of complex joystick movements or having to stop to use your seat lift? With our patented assignable buttons, you can assign C-Me to any button to lift or descend in one simple action. It's so easy you can even do it whilst driving!



Scan the QR code with your phone's camera to watch video.

STYLE

MAKE A STATEMENT

Treat every street like a catwalk with the Q300 R and its carefully selected range of 5 shroud colours. You can even compliment or contrast your choice with 7 colour accents. Scan the QR code below to find your perfect combination.



BLUE

PURPLE

YELLOW

RED

BLACK

built4me

Want something unique? Our Built-4-Me team can help you create a truly bespoke Q300 R that's out of this world! Ask your dealership for more information.

ACCENT COLOURS



BLACK



BLUE



GREEN



ORANGE



RED



SILVER



WHITE

TECHNICAL DATA

	Q300 R	Q300 R Mini
Seat width	385 - 520 mm	385 - 520 mm
Seat depth	410 - 510 mm	410 - 510 mm
Seat height	410 - 490 mm	410 - 490 mm
Back height	460 - 510 mm	460 - 510 mm
Lower leg length	305 - 460 mm	305 - 460 mm
Power tilt	0 -30°	0 -30°
Total width	580 mm	550 mm
Total length	From 1100 mm	From 1060 mm
Turning circle	1800 mm	1700 mm
Maximum safe slope	6°	6°
Maximum kerb climbing	100 mm	100 mm
Maximum range (ISO 7176-4)	32 km	32 km
Product weight (with 50 Ah batteries)	124 kg	124 kg
Maximum user weight	136 kg	136 kg
Drive wheels	14"	12"
Castor wheels	8" / 9"	8" / 9"
Maximum speed	10 kph	10 kph
Battery size	50 / 55 / 63 Ah	50 / 55 / 63 Ah
Electronics	VR2 / R-Net	VR2 / R-Net
Crash tested	Yes - approved for forward facing transport in a vehicle (ISO 7176-19). Attachment: 4 chair fixing points with 2 double straps rear and 2 single straps front. No anchored belt restraints allowed. The ease of access and manoeuvrability in vehicles depends on the wheelchair size.	
Transport without user	Fix the wheelchair safe in the vehicle (e.g. car or plane).	
Intended use and environment	For people with limited mobility, indoor and outdoor use. Wheelchair type class B	

Now that you've read all about the QUICKIE Q300 R,
why not arrange a free demonstration and experience
it for yourself?



Scan the QR code with your phone's camera
and book your free test-drive today.

Sunrise Medical cares about your safety. Please read your owner's manual for safe operating instructions and maintenance guidelines.

In line with our policy of continual product improvement, Sunrise Medical reserves the right to make changes at any time, and without notice. For the latest information, please visit www.sunrisemedical.co.uk

If you are visually impaired, this document can be viewed in PDF format at www.sunrisemedical.co.uk



Sunrise Medical Limited
Thorns Road, Brierley Hill, West Midlands, DY5 2LD
Tel.: +44 (0) 845 605 66 88
Fax: +44 (0) 845 605 66 89
E-Mail: enquires@sunmed.co.uk
www.sunrisemedical.co.uk