



REHASENSE



ICON 60

THE ACTIVE WHEELCHAIR



+44 (0) 333 220 0704



ukadmin@rehasense.com



www.rehasense.co.uk

ICON 60

THE ACTIVE WHEELCHAIR



Horizontal hand brake



Center of gravity adjustment

TECHNICAL DETAILS

Weight from	Transport weight	SWL (Safe Weight Limit)	Seat widths	Seat depths	Front frame angle (in deg) (to the seat)	Seat angles* (in deg)	Range of backrest angle (in deg)	Front seat heights (+ 3-5 cm cushion)	Rear seat heights (+ 3-5 cm cushion)	Footrest height (1.5 cm increment)	Center of gravity adjustment (1.3 cm increment)	Backrest height (6 cm of adjustment)	Rear wheels (2° camber & 1/2" axle)
9.0 kg	5.8 kg	120 kg	36-48 cm	42-46.5 cm (1.5 cm increment)	90°	6° 10° 13°	front: +20° back: up to -10°	47.5-50 cm (+ 3-5 cm cushion)	37.5-45 cm (+ 3-5 cm cushion)	4 positions (1.5 cm increment)	8 positions (1.3 cm increment)	6 cm of adjustment	2° camber & 1/2" axle

* depending on set up

DESIGNED WITH HIGH ADJUSTABILITY

The new ICON 60 manual wheelchair is one of the most adjustable wheelchairs currently available. This latest addition to the REHASENSE portfolio which features a lightweight rigid 7000 Series Aluminium Alloy frame and carbon fibre side guard has been designed to enable the user to adjust various elements of the chair to suit their individual needs. The frame length and seat angle are adjustable and for maximum comfort, the backrest - which is foldable to assist with transporting the ICON 60 by car - can be adjusted forward by up to 20° or backwards by up to 10°.

FEATURES AT A GLANCE


- Lightweight rigid frame wheelchair
- Seat width 36-50 cm, front width 26-40 cm, footrest width 19-33 cm
- Adjustable seat height (47,5-50 cm) and centre of gravity
- Front wheels angle adjustable
- Easily adjustable support bands for seat and backrest cushions
- Optional adjustable anti-tip

COLOUR

Black matt



SPECIAL FEATURES

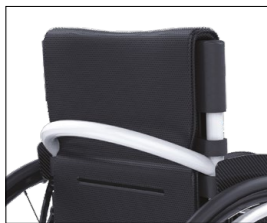

Tension and height adjustable backrest made with 3D breathable fabric



4" alu-core castor wheel with angle adjustable front fork



Lightweight carbon fibre side panel: foldable, removable and swing away




Height adjustable D style footplate



Low profile scissor lock

STANDARD CONFIGURATIONS

Item number	Product	Seat width	Seat height	Seat depth	Rear wheels	Camber	Axle	Colour
RW060036BK1	ICON 60	36 cm	49-51 cm	42-46.5 cm	24"	2°	1/2"	Black matt
RW060038BK1	ICON 60	38 cm	49-51 cm	42-46.5 cm	24"	2°	1/2"	Black matt
RW060040BK1	ICON 60	40 cm	49-51 cm	42-46.5 cm	24"	2°	1/2"	Black matt
RW060042BK1	ICON 60	42 cm	49-51 cm	42-46.5 cm	24"	2°	1/2"	Black matt
RW060044BK1	ICON 60	44 cm	49-51 cm	42-46.5 cm	24"	2°	1/2"	Black matt
RW060046BK1	ICON 60	46 cm	49-51 cm	42-46.5 cm	24"	2°	1/2"	Black matt
RW060048BK1	ICON 60	48 cm	49-51 cm	42-46.5 cm	24"	2°	1/2"	Black matt
coming soon	ICON 60	50 cm	49-51 cm	42-46.5 cm	24"	2°	1/2"	Black matt

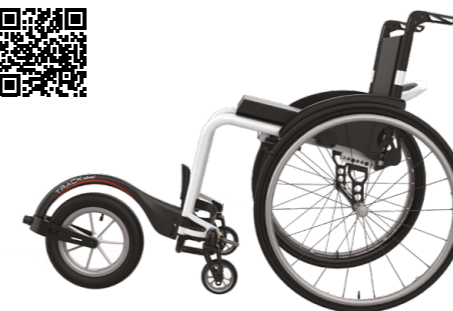

ICON 60 with TRACK WHEEL



ICON 60 with PAWS



REHASENSE

Rehasense UK Limited / Ireland
Unit 2, Harmony House
Treleigh Industrial Estate
Redruth, Cornwall TR16 4DE
United Kingdom



+44 (0) 333 220 0704



ukadmin@rehasense.com



www.rehasense.com

Follow us on social media



Rehasense's quality management complies with both the ISO 9001 and ISO 13485 standards, as confirmed by TÜV NORD. Rehasense's medical devices are CE-marked and meet the requirements of EU Regulation 2017/745 (MDR).

Copyright Rehasense© 2022



VI-2022