


Robby Bathing Aid

Simple and practical. Thanks to the low centre of gravity, the Robby Bathing Support needs comparatively little water in the bath. It can be folded up to a compact size allowing it to be stored easily. The compact size also makes it easy to take with you where needed.

Combining design with function, its four sets of ratchet joints and multiple positioning accessories support a wide range of sitting and lying positions for a safe and fun bath time.

What is special about me?

- Secure support while bathing
- Low body centre of gravity for optimum use of the bath
- Various sitting and lying positions
- Cover adaptable to the body shape with hook and loop straps
- Seat surface and back length grow with the child thanks to telescoping tubes and continuously adjustable supports

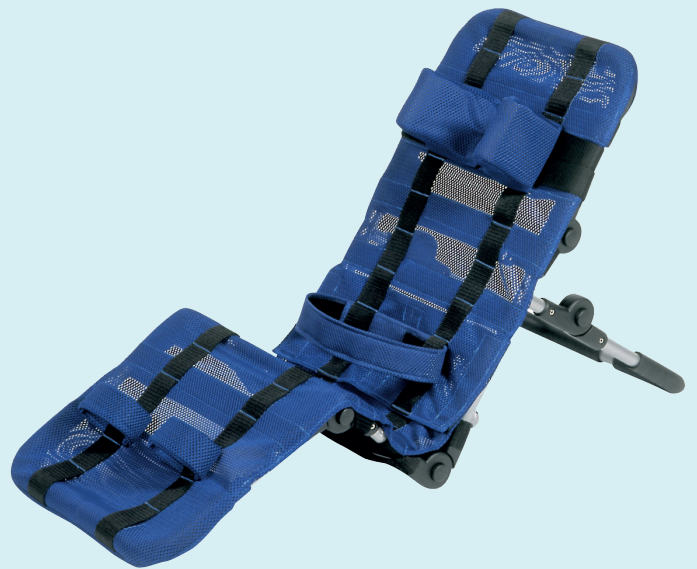

Robby Bathing Aid

When and where and I used?


- Bathing at home
- By the therapy or swimming pool
- Out and about at a family picnic or in the garden for temporary positioning



Adjustable straps for comfort



Easy folding for storage

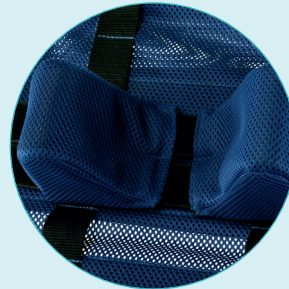

Headrest option



Other accessories available

Technical Data

	Size 1	Size 2
Max user weight	30kg	60kg
Seat depth (depending on strap tension and expansion of the telescoping tubes)	20 – 24cm	31 – 43cm
Backrest height	63 – 72cm	89 – 94cm
Total width	38cm	38cm
Total length (support surface)	115 – 133cm	147 – 173cm
Total weight (approx.)	3.8kg	4.3kg