


Easy to Fold




Step 1, please press the folding switch with one hand.

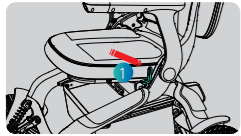


Step 2, push and lift the handle with the other hand until the Electric Wheelchair is fully folded.

Easy to unfold

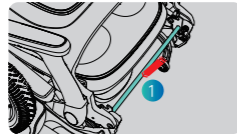

Step 1, press the folding bar down with one hand.

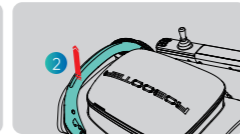

Step 2, pull the handle up with the other hand until the Electric Wheelchair is fully deployed.

PRODUCT PARAMETERS

Specification	E60 Pro	E60 Pro-A
Dimensions unfolded(L*W*H)	41.4"x25.2"x37.6"	41.4"x25.2"x37.6"
Dimensions folded(L*W*H)	33.7"x25.2"x20.7"	33.7"x25.2"x20.7"
Packaging Dimensions	36.6"x28.7"x21.9"	36.6"x28.7"x21.9"
Controller Height	30.3"	30.3"
Distance between armrests	18.9"	17.9-20.3" (stepless adjustable)
Seat Height	21"	21"
Seat Angle	4.5°	4.5°
Seat Width	16.5"	16.5"
Seat Depth	16.9"	16.9"
Backrest Height	18.1"	18.1"
Backrest Angle	18°	8-38°(4 adjustable angles of 10 ° per gear)
Footrest Height	6.1"	6.1"
Footrest Angle	5.5°	5.5°
Motor	2x150 W(Brushless motor)	2x150 W(Brushless motor)
Battery (Lithium-ion) 24V	20/25 Ah	20/25 Ah
Battery weight(20Ah)	7.3 lbs.	7.3 lbs.
Suspensions	Front wheel suspension	Front wheel suspension
Front Wheels Size (solid tire)	12 inch-Omni Wheels	12 inch-Omni Wheels
Rear Wheels Size (pneumatic tire)	12 inch	12 inch
Default Colors	Green\White\Grey	Green\White\Grey
App Control	Yes,remote control	Yes,remote control
Total mass	78.5 lbs.	82.5 lbs.
Max. User weight	330 lbs.	330 lbs.
Charging Time	7 h(20 Ah)、9 h(25 Ah)	7 h(20 Ah)、9 h(25 Ah)
Range (up to)	14.3 miles(20 Ah)、18 miles(25 Ah)	14.3 miles(20 Ah)、18 miles(25 Ah)
Options-Maximum Speed	5 mph,6.2 mph	5 mph,6.2 mph
Obstacle Climbing	3.1"	3.1"
Turning Radius	< 33.5"	< 33.5"
Static Stability (Climbing angle)	15°	15°
Waterproof Level	IPX5	IPX5

*The above data is obtained by manual measurement and may have some errors. Please refer to the actual product;

*Depending on production date of electric wheelchair, product hardware and software may vary slightly, the specific situation and functions will be subject to real world situations and performance.

REBOOT YOUR LIFE

COLOR



ROBOOTER

REBOOT YOUR LIFE

E60 Pro /E60 Pro-A

Off-road power wheelchair



OUR FULL PRODUCT LINE


■ Standing Wheelchairs for Rehabilitation



XZ-Droid Smart




X40

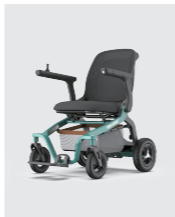

E40




E60-A

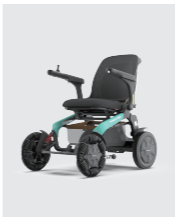

E60 Pro-A

■ Mobility aids ■

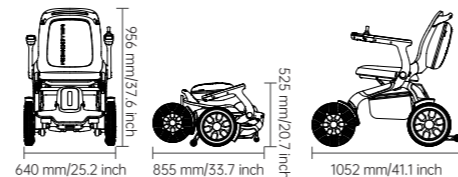

Shanghai Bangbang Robotics Co.,Ltd
Floor 5, Building 3, No.188 Zhongchen Road,
Songjiang District, Shanghai, China.
www.robooter.com

E60 Pro/E60 Pro-A

E60 Pro/E60 Pro-A is the first of its kind, smart, off-road power wheelchair launched by ROBOOTER. It delivers high endurance and fast driving speed, we think you will agree that it is the perfect blend of design and function.



Magnesium Aluminum Alloy

Light weight, loading 150kg.



Multifunctional armrest

Built-in front lighting, USB/USBC charging port.




Storage designed to go with you

Water-repellent large-capacity backpack and mesh pocket, loading 10kg.

BUILT TO PROVIDE THE ULTIMATE IN OFF-ROAD DRIVING

Front independent suspension system

Adjustable suspension spring with strong spring coils, ride in comfort.

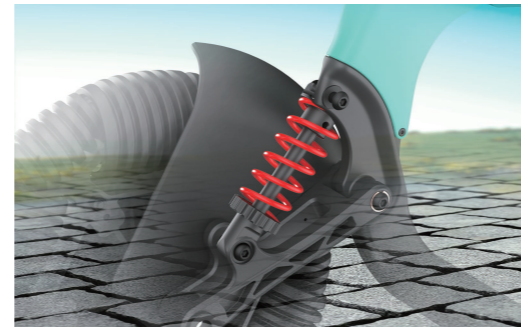

12-inch front and rear wheels, obstacle clearance height up to 3.1-inch.

The front wheels are omnidirectional and can move laterally, making it easy to cross and turn in narrow spaces.



3.1-inch
Maximum obstacle crossing height

5.9-inch
Maximum ditch width

15°
Climbing angle

33.5-inch
Minimum turning radius


Adjustable back

The E60 Pro-A back angle can reach 8-38°, relieving back fatigue and soreness.



Adjustable seat width

Adjust the width between the armrests from 17.9-20.3", for personalize comfort.

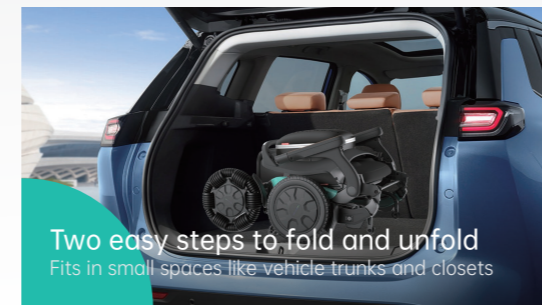


Two easy steps to fold and unfold
Fits in small spaces like vehicle trunks and closets



App with Bluetooth Connectivity
Adjust setting, monitor device, and even remotely drive!



BMS intelligent management system
for lithium batteries

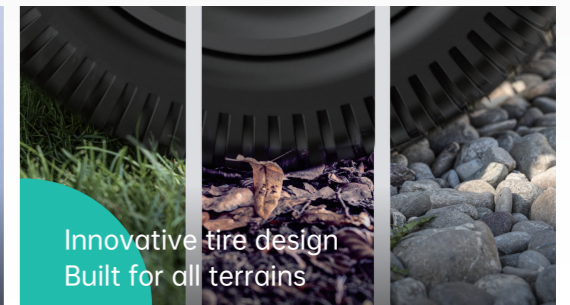

Innovative tire design
Built for all terrains