

 **loopwheels** **urban**

*It's the only wheel you'll need*

Go light  
Go comfortable  
Go anywhere

[loopwheels.com](http://loopwheels.com)





# Lightweight, vibration-reducing wheels for everyday use.

Safe levels of whole body vibration for over 2 hours every day.

## ISO SAFETY MEASURES 2631-1 1985

This measures the exposure to vibration in a seated position, before action is needed to protect health.


### ATV (average daily use) tested at 12km/h over uneven surfaces

Spinergy 24" day chair wheel	safe for 1 hour 3 minutes
Softwheel 24" day chair wheel	safe for 1 hour 18 minutes (+24%)
<b>Loopwheels Urban</b>	<b>safe for 2 hours 10 minutes (+106%)</b>

Loopwheels Urban is safe for over **TWICE THE TIME** vs a Spinergy wheel on off-road surfaces at an average of 12km/h.

## PRODUCT SPECIFICATION

- Wheel size:** 24" (540mm)
- ETRTO number:** 25-540 (24x1.00)
- Axle bearing:** 1/2" (12.7mm) and 12mm
- Push rim:** 16 rivnut fitting
- Rivnut centres measurement:** 515mm
- Spring/loop materials:** Carbon fibre composite
- Rim:** Double-walled aluminium, black
- Inner rim width:** 16mm
- Hub:** CWC machined aluminium, anodised black sealed cartridge bearings
- Bearing width at faces:** 1.8" (45.8mm)
- Wheel weight:** 1.3kg (2.87lb)  
NB Wheel weights without push rim, tyre, or axle
- Maximum load (user and chair):** 120kg (265lb)
- Recommended tyres:** Schwalbe Marathon Plus

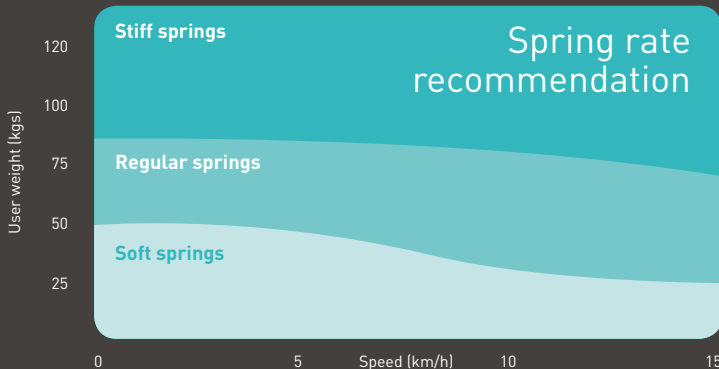

Loopwheels Urban is available in 3 spring strengths – **Soft springs**, **Regular springs** and **Stiff springs**. Use the chart above to estimate the most usual higher speed on the wheels and the fully-loaded weight of the chair.

## Jelly Products Ltd

Unit 202 Boughton Industrial Estate North,  
Newark, Nottinghamshire NG22 9LD  
T: +44 (0) 1623 860909  
E: info@loopwheels.com  
W: loopwheels.com

To learn how Loopwheels can make the difference, visit [loopwheels.com](https://loopwheels.com)